overfibers
improve your performance

Made in Italy
by Overfibers S.r.l


CE 0425

GB	Glass fiber endodontic posts	Instructions for use
I	Perni endocanalari in fibra di vetro	Istruzioni d'uso
F	Tenons endocanalaires en fibre de verre	Instructions d'emploi
E	Postes endo-radiculares de fibra de vidrio	Instrucciones de uso
PL	Wkłady korzeniowe z włókna szklanego	Instrukcja użytkowania


Duo Over Post

Cementation technique


Overfibers Post clinical Use

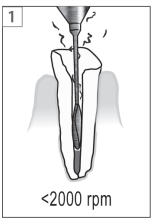
Canal Disinfection


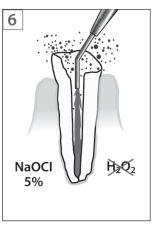
A. PROCEDURE FOR OVERCEM SA AND ALL SELF-ADHESIVE CEMENTS

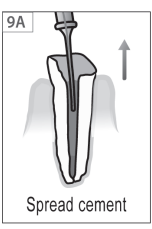B. PROCEDURE FOR RESIN CEMENTS & ADHESIVES

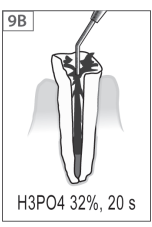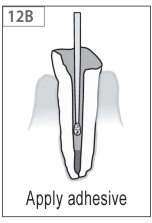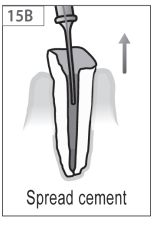Overfibers S.r.l. Via Malatesta, 7 - 40026 Imola (BO) - Italy
Tel. +39 0542.52153 - Fax +39 0542.1900078
e-mail: info@overfibers.com www.overfibers.com

Duo Over Post

English

Translucent and radiopaque glass fiber endodontic posts

Duo Over Posts with an extra thin point were designed to be used in pairs in "ribbon" canals in which the width, for example the vestibule-lingual one, strongly dominates the mesio-distal one. They are composed of high-resistance fiberglass made of an epoxy resin matrix. The translucence of the post material allows the realization of highly aesthetic restorations. Their elastic modulus and durability characteristics are calculated so as to reduce the risk of fractures at the root and at the same time ensure the maximum stability of the restoration. The high surface roughness maximizes the adhesion of the cement and the composite material.

Sizes and color codes

Duo Over Posts are available in two sizes of 19 mm in length, each one identified by a color code:

Duo Over Post			
Color	Ø Max	Ø Min	Taper (%)
Yellow #1	1.5	0.50	0.071
Red #2	1.8	0.60	0.086

Tab. I

Indications
Duo Over Posts have been designed to reconstruct teeth with flared or ribboned channels, in particular premolars, canines, or superior central incisors. They should be used in pairs, and if the clinical case so requires, even in threes. The two sizes can be combined with each other in the same canal.

Steps for clinical use
1. The post should be used after having successfully completed treatment of the canal and having removed any trace of dentin decay;
2. To calculate the post's insertion depth, follow the X-rays or data acquired during canal treatment;
3. Completely remove the material sealing the canal until the calculated depth. Minimize removal of canal dentin. At least 3 - 4mm of the apical seal must remain to guarantee endodontic success.
4. For prosthetic restorations: at least 2 mm of dentin must remain above the foreseen margin of coronal closure to guarantee the "ferule" effect;
5. Select the post to be used based on the type of tooth and the width of the ribbon-like canal;
• the post must be the largest possible compatible with the canal diameter;
• the portion that enters the canal must be at least as long as that which remains outside upon restoration;
6. Irrigate the canal with detergent or disinfectant solutions (e.g. 5% NaOCl, EDTA);
7. Application procedures: Check the post's fit in the canal. The post must reach the previously determined depth and have minimal play once inserted (2-3 tenths of a mm). If the play is greater, try the next size post.
• If necessary, shorten the coronal portion of the post, bringing it to the desired length. Use only diamond rotating instruments of medium or large granulometry, and do not use cutting instruments such as scissors or side cutters so as not to damage the material. Carry out this operation far from the patient, using gloves, a mask, and protective goggles.
• Cleanse the post by immersing it in pure ethyl alcohol or coating it with 32% phosphoric acid gel for 60 seconds; then wash it with plenty of water and dry it.
• As the Overfibers post surfaces are particularly rough, salinization treatments are optional.
• We recommend using a resin cement as OverCEM SA for the

Duo Over Post

English

cementing of Over Posts to insure the best possible adhesion to the root canal.
• If a dual system is chosen, the coronal extremity of the post must be wet with adhesive so as to increase translucency and therefore the passage of light from the photo-polymerizing lamp. Keep the light source very close and aligned with the long axis of the post for at least 60 seconds.
• Reconstruction of the stump: choose the desired reconstructive system; the dentin adhesive treatment must be conducted according to the manufacturer's instructions. The use of a reconstruction composite specific for stumps is recommended. Place the composite directly on the external part of the post, completely enveloping it, and being careful not to create empty spaces. Follow the manufacturer's instructions in this case as well.
• CAUTION: temporary cements or disinfectant substances containing phenols (e.g. eugenol) previously applied to a tooth can greatly reduce the bond to the dentin and the strength of resinous cements and reconstruction composites!
• To create the ferule effect, ensure that the preparation margin is set to at least 2mm in healthy tissue; on teeth prepared for the anchorage of removable prostheses (RPD), such a margin should be greater if possible.
• After the preparation, apply the temporary, avoiding cements or disinfectants that contain phenols (e.g. ZOE cements).
8. Removal of the post in the case of reintervention (e.g. restoration fracture, secondary cavities, canal retreatment): use diamond rotating or bladed instruments with extreme caution to avoid perforations. The special removal tool kit available commercially can be used. The use of magnifying lenses is recommended.

Precautions
• To obtain the safety data sheets (Overfibers MSDS), go to the www.overfibers.com site

Storage
The product must be conserved in a cool, dry place protected from sunlight.
The device is disposable.

Disposal
Dispose in accordance with the regulations in force.

Warnings and contraindications
Do not use in patients with known or suspected allergies against the materials constituting the device.
Duo Over Post
• The posts can withstand a sterilization procedure in envelopes through autoclaving at 135°C and 2.1 bar of pressure for 5 minutes. This procedure must be applied in all cases of suspected or established biological contamination of the post.
• Avoid the prolonged and repeated breathing of dust arising from the abrasion of the device during clinical procedures. The use of protective masks, high-velocity suction, and rubber dams is recommended.
• Keep away from children.
• The product is intended solely for dental use.

Manufacturer: Overfibers S.r.l
Via Malatesta, 7 - 40026 - Imola (Bo) Italy

1) This material has been developed solely for use in dentistry. It must be employed by exclusively following the instructions for use. The manufacturer assumes no liability for damages resulting from the failure to observe the instructions for use or for a use other than the scope foreseen for this product. Users are responsible for testing the material in usage fields not explicitly mentioned in the instructions for use. The descriptions and data do not constitute any guarantee of the product's properties and are not binding.

Production of the instructions for use: 03/2021, Rev. 09

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Producent: Overfibers S.r.l - Via Malatesta, 7 - 40026 - Imola (Bo) Italy