overfibers
improve your performance

Made in Italy
by Overfibers S.r.l


CE 0425


GB	Glass fiber endodontic posts	Instructions for use
I	Perni endocanalari in fibra di vetro	Istruzioni d'uso
F	Tenons endocanalaires en fibre de verre	Instructions d'emploi
E	Postes endo-radiculares de fibra de vidrio	Instrucciones de uso
PL	Wkłady korzeniowe z włókna szklanego	Instrukcja użytkowania


Prosthetic Over Post


Cementation technique


Overfibers Post clinical Use

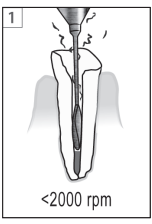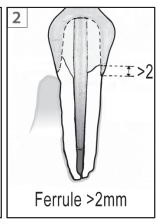
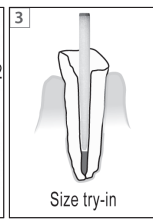
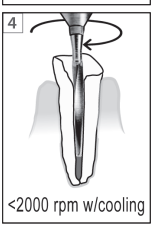
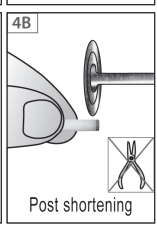
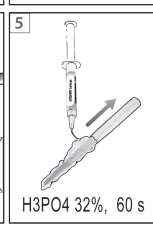
Canal Disinfection

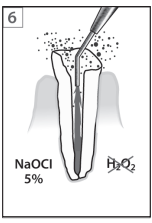
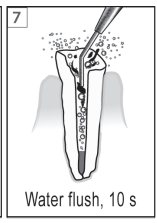
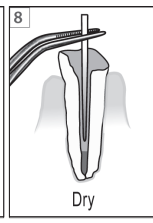
A. PROCEDURE FOR OVERCEM SA AND ALL SELF-ADHESIVE CEMENTS


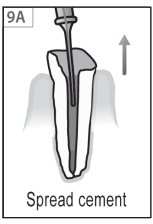
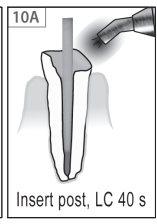
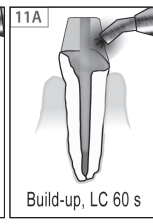
B. PROCEDURE FOR RESIN CEMENTS & ADHESIVES

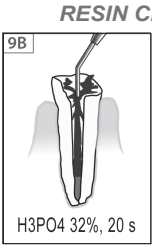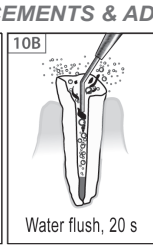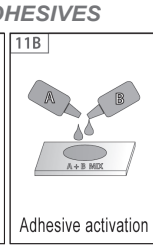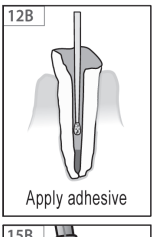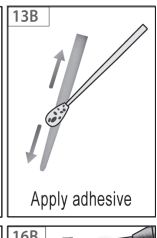
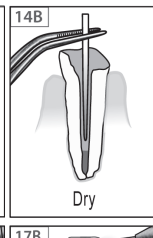
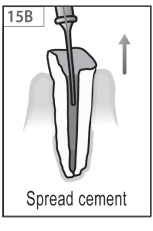
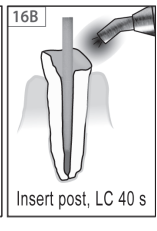
Overfibers S.r.l. Via Malatesta, 7 - 40026 Imola (BO) - Italy
Tel. +39 0542.52153 - Fax +39 0542.1900078
e-mail: info@overfibers.com www.overfibers.com

Prosthetic Over Post
Translucent and radiopaque glass fiber endodontic posts.
The conical-cylindrical shaped Prosthetic Over Posts are particularly well suited to teeth with little or no residual dentin for treatment with prosthetic reconstructions (crowns or bridges), as they achieve the maximum diameter of emergence from the root. They are composed of high-resistance fiberglass in an epoxy resin matrix. The translucency of the post material allows the realization of highly aesthetic restorations. Their elastic modulus and durability characteristics are calculated so as to reduce the risk of fractures at the root and at the same time ensure the maximum stability of the restoration. The high surface roughness maximizes the adhesion of the cement and the composite material.

Sizes and color codes

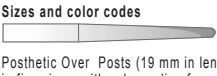Posthetic Over Posts (19 mm in length) and Prosthetic Drills, are available in five sizes, with color coding for a corresponding use:

Prosthetic Over Post				Prosthetic Drill
Color	Max Ø	Min Ø	Taper (%)	Color
Black #0	1,20	0,70	0,04	Black #0
Yellow #1	1,40	0,80	0,05	Yellow #1
Red #2	1,60	0,80	0,066	Red #2
Blue #3	1,80	0,90	0,08	Blue #3
Green #4	2,10	1,10	0,08	Green #4

Tab. I

Indications
Prosthetic Over Posts have been designed to suit all clinical cases, and in particular teeth that are extensively destroyed and very stressed. The available sizes allow all dental typologies to be treated (Tab. II).

Post Code	Teeth
Black #0	Inferior incisors
Yellow #1	Inferior and superior incisors, first superior premolars
Red #2	Superior incisors, canines, premolars, molars
Blue #3	Superior incisors, canines, premolars, molars
Green # 4	Superior incisors, canines, premolars, molars

Tab. II: General indications

1.The post should be used after having successfully completed treatment of the canal and having removed any trace of dentin decay;
2.To calculate the post's insertion depth, follow the X-rays or data acquired during canal treatment;
3.Completely remove the material sealing the canal until the calculated depth using endodontic drills such as Gates-Glidden® or Largo® at low velocity (1500-2000 RPM). Minimize removal of canal dentin. At least 3 - 4mm of the apical seal must remain to guarantee endodontic success.
4. For prosthetic restorations: at least 2 mm of dentin must remain above the foreseen margin of coronal closure to guarantee the "ferule" effect;
5. Select the post to be used based on the type of tooth and the width of the canal (Tab. II);
• the post must be the largest possible compatible with the canal diameter;
• the portion that enters the canal must be at least as long as that which remains outside upon restoration;
6. Once the post has been chosen, prepare the canal with a suitable cutter Prosthetic Drill, refer to table I.
Prosthetic Drill must be sterilized prior to use.
7. Irrigate the canal with detergent or disinfectant solutions (e.g. 5% NaOCl, EDTA);
8. Application procedures: Check the post's fit in the canal. The post must reach the previously determined depth and have minimal play once inserted (2-3 tenths of a mm). If the play is greater, try the next size post.
• If necessary, shorten the coronal portion of the post, bringing it to the desired length. Use only diamond rotating instruments of medium or large

Prosthetic Over Post English

granulometry, and do not use cutting instruments such as scissors or side cutters so as not to damage the material. Carry out this operation far from the patient, using gloves, a mask, and protective goggles.
• Cleanse the post by immersing it in pure ethyl alcohol or coating it with 32% phosphoric acid gel for 60 seconds; then wash it with plenty of water and dry it.
• As the Overfibers post surfaces are particularly rough, silanization treatments are optional.
• We recommend using a resin cement as OverCEM SA for the cementing of Over Posts to insure the best possible adhesion to the root canal.
• If a dual system is chosen, the coronal extremity of the post must be wet with adhesive so as to increase translucency and therefore the passage of light from the photo-polymerizing lamp. Keep the light source very close and aligned with the long axis of the post for at least 60 seconds.
• Reconstruction of the stump: choose the desired reconstructive system; the dentin adhesive treatment must be conducted according to the manufacturer's instructions. The use of a reconstruction composite specific for stumps is recommended. Place the composite directly on the external part of the post, completely enveloping it, and being careful not to create empty spaces. Follow the manufacturer's instructions in this case as well.
• CAUTION: temporary cements or disinfectant substances containing phenols (e.g. eugenol) previously applied to a tooth can greatly reduce the bond to the dentin and the strength of resinous cements and reconstruction composites!
• To create the ferule effect, ensure that the preparation margin is set to at least 2mm in healthy tissue; on teeth prepared for the anchorage of removable prostheses (RPD), such a margin should be greater if possible.
• After the preparation, apply the temporary, avoiding cements or disinfectants that contain phenols (e.g. ZOE cements).
9. Removal of the post in the case of reintervention (e.g. restoration fracture, secondary cavities, canal retreatment): use diamond rotating or bladed instruments with extreme caution to avoid perforations. The special removal tool kit available commercially can be used. The use of magnifying lenses is recommended.

Precautions
• To obtain the safety data sheets (Overfibers MSDS), go to the www.overfibers.com site

Storage
The product must be conserved in a cool, dry place protected from sunlight. The device is disposable.

Disposal
Dispose in accordance with the regulations in force.

Warnings and contraindications
Do not use in patients with known or suspected allergies against the materials constituting the device.
• Prosthetic Over Post can withstand a sterilization procedure in envelopes through autoclaving at 135°C and 2.1 bar of pressure for 5 minutes. This procedure must be applied in all cases of suspected or established biological contamination of the post.
• Avoid the prolonged and repeated breathing of dust arising from the abrasion of the device during clinical procedures. The use of protective masks, high-velocity suction, and rubber dams is recommended.
• Keep away from children.
• The product is intended solely for dental use.

Manufacturer: Overfibers S.r.l - Via Malatesta, 7 - 40026 - Imola (Bo) Italy
1) This material has been developed solely for use in dentistry. Use of this material must adhere strictly to instructions. The manufacturer shall assume no liability for damage resulting from failure to observe the instructions for use or for any use other than that intended for this product. Users are responsible for testing the material in any usage fields not explicitly mentioned in the instructions. None of the descriptions or data constitute any guarantee of product properties and are not binding.

Production of the instructions for use: 01/2021, Rev. 09

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Producent: Overfibers S.r.l - Via Malatesta, 7 - 40026 - Imola (Bo) Italy